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## Instructions for Use

189826

EN	Instructions for Use	2
SL	Navodila za uporabo	11
BG	Инструкции за употреба	20
HR	Upute za uporabu	29
SK	Návod na použitie	38
HU	Használati útmutató	47
EL	Οδηγίες χρήσης	56
LV	Lietošanas pamācība	65
LT	Naudojimo instrukcija	74
ET	Kasutusjuhend	83

# Contents

Contents .....	2
1 Description and Intended Purpose .....	3
2 Safety Information .....	4
3 Construction .....	5
4 Function.....	6
5 Maintenance .....	6
6 Limitations on Use .....	7
7 Fitting Advice.....	7
8 Assembly Instructions .....	8
9 Technical Data .....	9
10 Ordering Information .....	9

---

# 1 Description and Intended Purpose

These Instructions for Use are intended for use by the practitioner and user unless otherwise stated.

The term *device* is used throughout these instructions for use to refer to M30 Mini BladeXT Adapter.

Please read and ensure you understand all instructions for use, in particular all safety information and maintenance instructions.

## Application

This device is to be used exclusively as part of a lower limb prosthesis.

Intended for a single user.

The device allows connection of a Mini BladeXT to a socket fitted with a Blatchford M30 4-Prong Laminate Adapter 189626, via single bolt distally and via a threaded section to the laminate adapter.

## Feature

- Lightweight, strong, titanium construction

## Activity Level

This device is suitable for Activity Levels 1-4 (weight limits apply, see Section 9 *Technical Data*). Of course there are exceptions and in our recommendation we want to allow for unique, individual circumstances and any such decision should be made with sound and thorough justification.

### Activity Level 1

Has the ability or potential to use a prosthesis for transfers or ambulation on level surfaces at fixed cadence. Typical of the limited and unlimited ambulator.

### Activity Level 2

Has the ability or potential for ambulation with the ability to traverse low level environmental barriers such as curbs, stairs, or uneven surfaces. Typical of the limited community ambulator.

### Activity Level 3

Has the ability or potential for ambulation with variable cadence.

Typical of the community ambulator who has the ability to traverse most environmental barriers and may have vocational, therapeutic, or exercise activity that demands prosthetic utilization beyond simple locomotion.

### Activity Level 4

Has the ability or potential for prosthetic ambulation that exceeds basic ambulation skills, exhibiting high impact, stress, or energy levels. Typical of the prosthetic demands of the child, active adult, or athlete.

## Contraindications

There are no known contraindications if this device is used in accordance with these instructions.

## Clinical Benefits

- Allows connection of the Mini BladeXT to socket

---

## 2 Safety Information

**This warning symbol highlights important safety information which must be followed carefully.**

Any changes in the performance or function of the limb, e.g. movement, should be immediately reported to your service provider.

Always use a hand rail when descending stairs and at any other time if available.

The device is not suitable for extreme sports, running or cycle racing, ice and snow sports, extreme slopes and steps. Any such activities undertaken are done so completely at the users' own risk. Recreational cycling is acceptable.

The device is designed for prolonged submersion and suitable for immersion in fresh water only. Ensure any use of the device in water complies with the conditions given in Section 6 *Limitations on Use*.

Be aware of finger trap hazard at all times.

Assembly, maintenance and repair of the device must only be carried out by a suitably qualified practitioner.

The user must not adjust or tamper with the setup of the device.


The user should be advised to contact their practitioner if their condition changes.

# 3 Construction

## Principal Parts

- M30 Adapter Titanium
- Bolt Titanium
- Alignment Washer Stainless steel
- Alignment Nut Stainless steel

## Component Identification

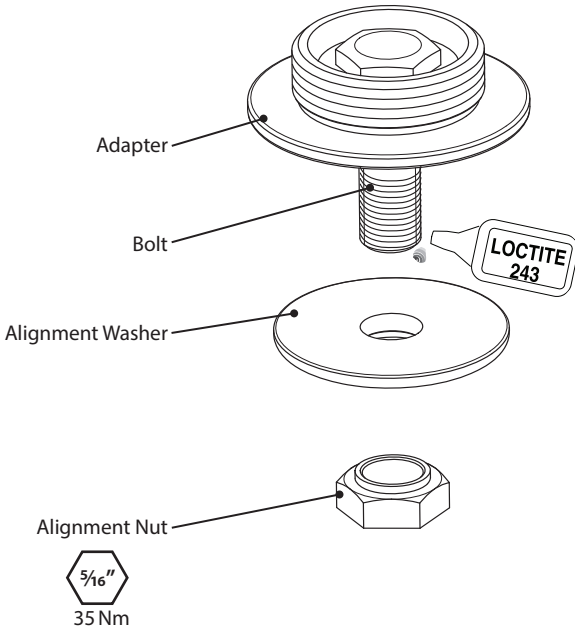---

## 4 Function

This device provides shift along the sagittal plane, which allows for a customized limb alignment to suit the user. The device also allows rotational adjustment at the same single bolt attachment, and offers alternative rotation adjustment possibilities when interfacing the threaded section of the adapter.

When the nut is tightened securely using Loctite 243 and the torque setting shown in the *Construction* section, this device offers a secure way to connect two parts of the limb.

---

## 5 Maintenance

Visually check the device regularly.

Report any changes in performance of this device to the practitioner/service provider e.g. unusual noises, loss of alignment or looseness.

Inform the practitioner/service provider of any changes in body weight and/or activity level.

### **Cleaning**

Use a damp cloth and mild soap to clean outside surfaces. **Do not** use aggressive cleansers.

### ***The remaining instructions in this section are for practitioner use only.***

This maintenance must be carried out only by competent personnel (practitioner or suitable trained technician).

The following routine maintenance is to be carried out at least annually:

- Check the alignment bolt for tightness, see *Construction* section. If loose, remove and clean, then apply Loctite and tighten to the correct torque setting.
- Check for defects that could affect proper function.

Ensure the user has read and understood all safety and user-level maintenance information.

Advise the user that a regular visual check of the device is recommended and signs of wear that may affect function should be reported to their service provider.

Advise the user to inform the practitioner/service provider of any changes in body weight and/or activity level.

If this device is used for extreme activity, the maintenance level and interval should be reviewed and if required advice and technical support sought to plan a new maintenance schedule dependent upon the frequency and nature of the activity. This should be determined by a local risk assessment carried out by a suitably qualified individual.

---

## 6 Limitations on Use

### Intended Life

A local risk assessment should be carried out based upon activity and usage.

### Lifting Loads

User weight and activity is governed by the stated limits.

Load carrying by the user should be based on a local risk assessment.

### Environment

This product is waterproof up to a maximum depth of 1 meter. To prevent wear or damage, thoroughly rinse with fresh water after use in abrasive environments such as those that contain sand or grit. Thoroughly rinse with fresh water after use in salt or chlorinated water.

Exclusively for use between -15 °C and 50 °C  
(5 °F to 122 °F).

Suitable for submersion

---

## 7 Fitting Advice


*The instructions in this section are for practitioner use only.*

Symptom	Solution
A recurring noise occurs.	The bolt must be tightened to the torque ratings specified in this document. Apply Loctite.
The device moves out of position.	Advise the user they must not use the device until adjusted, repaired or replaced. Inspect device. If components are worn or damaged, then remove and replace with new components. Reposition as necessary and retighten as per <i>Construction</i> section..

# 8 Assembly Instructions

Before you begin assembly, make sure that the interfacing surface of the MiniBladeXT is clean and free from grease (clean as necessary).

Assemble as per the diagram below to the torque setting shown. Apply Loctite 243 to the exposed threads before fitting the nut.

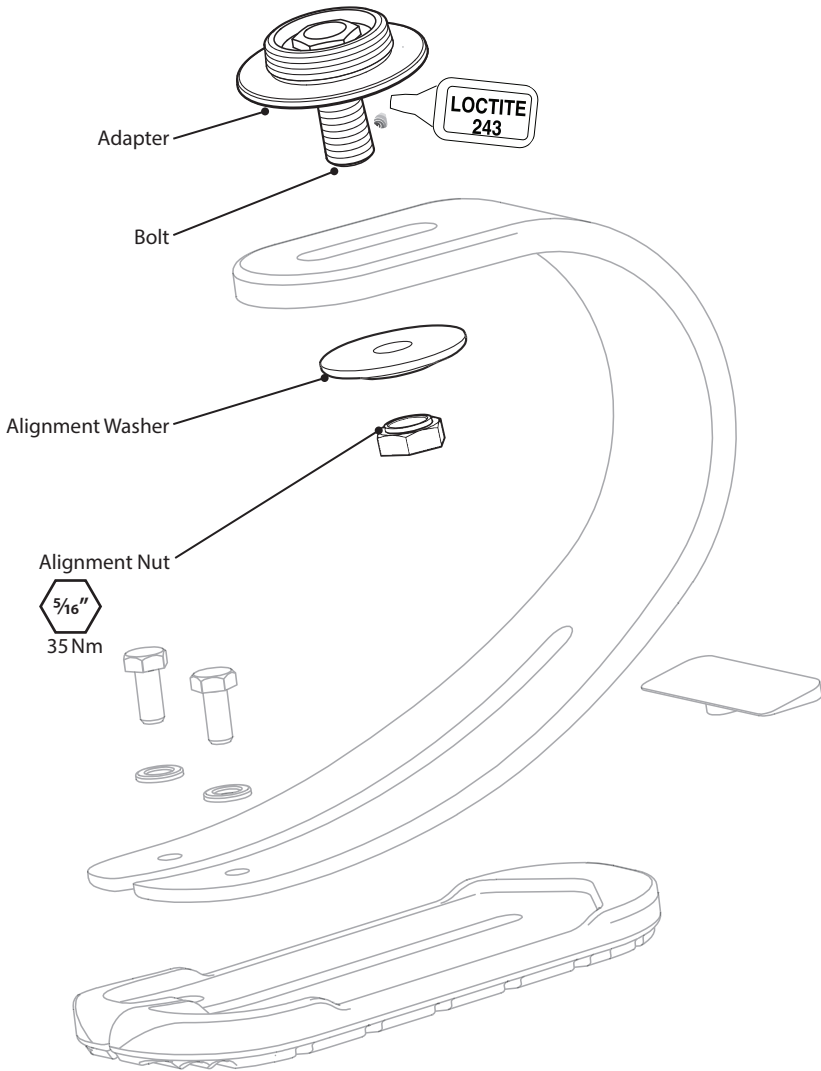
---

## 9 Technical Data

Operating and Storage Temperature Range:	-15 °C to 50 °C (5 °F to 122 °F)
Component Weight:	80 g (2.8 oz)
Activity Level:	1–4
Maximum User Weight:	60 kg (132 lb)
Range of Adjustment	360° rotation ±15 mm A-P shift
Build Height:	12 mm
Proximal Attachment:	M30 4 Prong Laminate Adapter
Distal Attachment:	Mini Blade XT

---

### Fitting Length
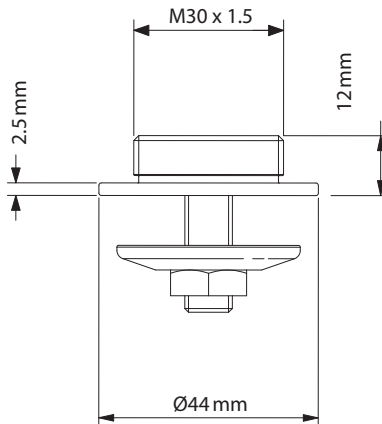
## 10 Ordering Information

**Order Number 189826**

---

## Liability

The manufacturer recommends using the device only under the specified conditions and for the intended purposes. The device must be maintained according to the instructions for use supplied with the device. The manufacturer is not liable for any adverse outcome caused by any component combinations that were not authorized by them.

## CE Conformity

This product meets the requirements of the European Regulation EU 2017/745 for medical devices. This product has been classified as a class I device according to the classification rules outlined in Annex VIII of the regulation. The EU declaration of conformity certificate is available at the following internet address: [www.blatchford.co.uk](http://www.blatchford.co.uk).

Medical Device

Single Patient – multiple use

## Compatibility

Combination with Blatchford branded products is approved based on testing in accordance with relevant standards and the MDR including structural test, dimensional compatibility and monitored field performance.

Combination with alternative CE marked products must be carried out in view of a documented local risk assessment carried out by a practitioner.

## Warranty

This device is warranted for 24 months.

The user should be aware that changes or modifications not expressly approved could void the warranty, operating licenses and exemptions.

See the Blatchford website for the current full warranty statement.

## Reporting of Serious Incidents

In the unlikely event of a serious incident occurring in relation to this device it should be reported to the manufacturer and your national competent authority.

## Environmental Aspects

This product is constructed from recyclable materials. Where possible the components should be recycled in accordance with local waste handling regulations.

## Retaining the Packaging Label

The practitioner is advised to keep the packaging label as a record of the device supplied.

## Trademark Acknowledgements

Blatchford is a registered trademark of Blatchford Products Limited.

## Manufacturer's Registered Address

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**Blatchford Europe GmbH**
**Endolite India Ltd.**